



Activator

Instructions for Use

Scan the QR



Intended use

The Activator is required for both the activation and deactivation of the **Miniaturized Electromagnetic Device** (MED).

Introduction

The Activator switches the MED from the sleep mode to the active mode:

Sleep mode– while in storage, the MED will not produce an electromagnetic field. You cannot install the MED in this mode.

Active mode– the MED will produce an electromagnetic field. You can only install the MED in this mode.


Indications

The Activator is required for activation of all MED devices.

Note

Read the instructions for the MED before using the Activator.

Caution

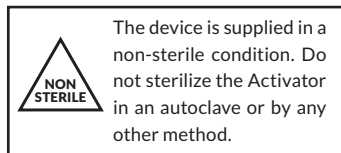

Manufacturing date



The product's year and month of manufacturing date appears on the product label located on the bottom of the Activator.

Point of Use

The health care institution is the point of use.

Parts and functions



1. The MED is placed in the activation socket for activation/deactivation while it is still in its inner vial.

2. The activation button is used to switch the MED between the 2 modes:

(A) Sleep mode:

the activation button will light up in red for 15 seconds to indicate that the MED is in the sleep mode.

(B) Active mode:

the activation button will light up in blue for 15 seconds to indicate that the MED is in the active mode.

3. Low battery: *the low battery indicator will blink in red to indicate that the batteries need to be replaced. The Activator uses two 1.5V AAA batteries.*

Packaging:

The package contains one ACTIVATOR and two 1.5V AAA batteries.





Activator

Instructions for Use

Scan the QR



First use

1. Open the Activator box, pull the tab and open the upper container.
2. Insert the two 1.5V AAA batteries (supplied with the Activator) into the battery housing on the bottom of the Activator.

Activation

STEP 1 - Place the inner vial containing the MED in to the activation socket. *The activation button will light up in red*

NOTE: The MED can only be taken out of the inner vial after activation and just before manually threading it into the implant.

STEP 2 - Press the activation button to begin the activation process. During this time, *the activation button will blink in blue* for 2 seconds while the Activator emits a *beeping sound*.

After the activation process is over, *the activation button will light up in blue* to indicate the MED is generating an electromagnetic field and is ready for use in patients.

Do not use the MED if the activation button does not *light up in blue* after the activation process. This may indicate that either the MED or the Activator may be defective. In this case, contact our customer service.

STEP 3 - Remove the inner vial that contains the MED from the activation socket and follow the MED instructions for use.

Deactivation

STEP 1 - Place the inner vial containing the MED in the activation socket. *The activation button will light up in blue.*

STEP 2 - Press the activation button for 3 seconds. During the deactivation process, *the activation button will blink in red* for 2 seconds, and the Activator will emit a *long beeping sound*.

After the deactivation process is over, *the activation button will light up in red* for 15 seconds to indicate that the MED is in the sleep mode.

First use

1. Store the device at room temperature.
2. Do not place the device in contact with water or any other liquids.
3. Handle with care.

Manufacturer's Information:



Magdent Ltd.

B.S.R Tower No. 3, 15th floor
9 Mesada st. Bnei-Brak, 5120109, Israel
Tel: +972-3-373-1701
Email: support@magdentmed.com

For electromagnetic safety and compliance information, please visit our website:

www.magdentmed.com

EC

REP

MedNet GmbH, Borkstrasse 10,
48163 Muenster, Germany



CE Marking



Double insulated
appliance



Follow instructions for
use on the website



Do not throw away -
Electrical equipment